

TECNO SAVE


 **Tecnoline**
Medical & Surgical Solutions

POST-OPERATIVE ORTHOPAEDIC BLOOD RECOVERY SYSTEM




ISO 13485

SISTEMA DI GESTIONE QUALITÀ UNI CEI EN ISO 13485:2016

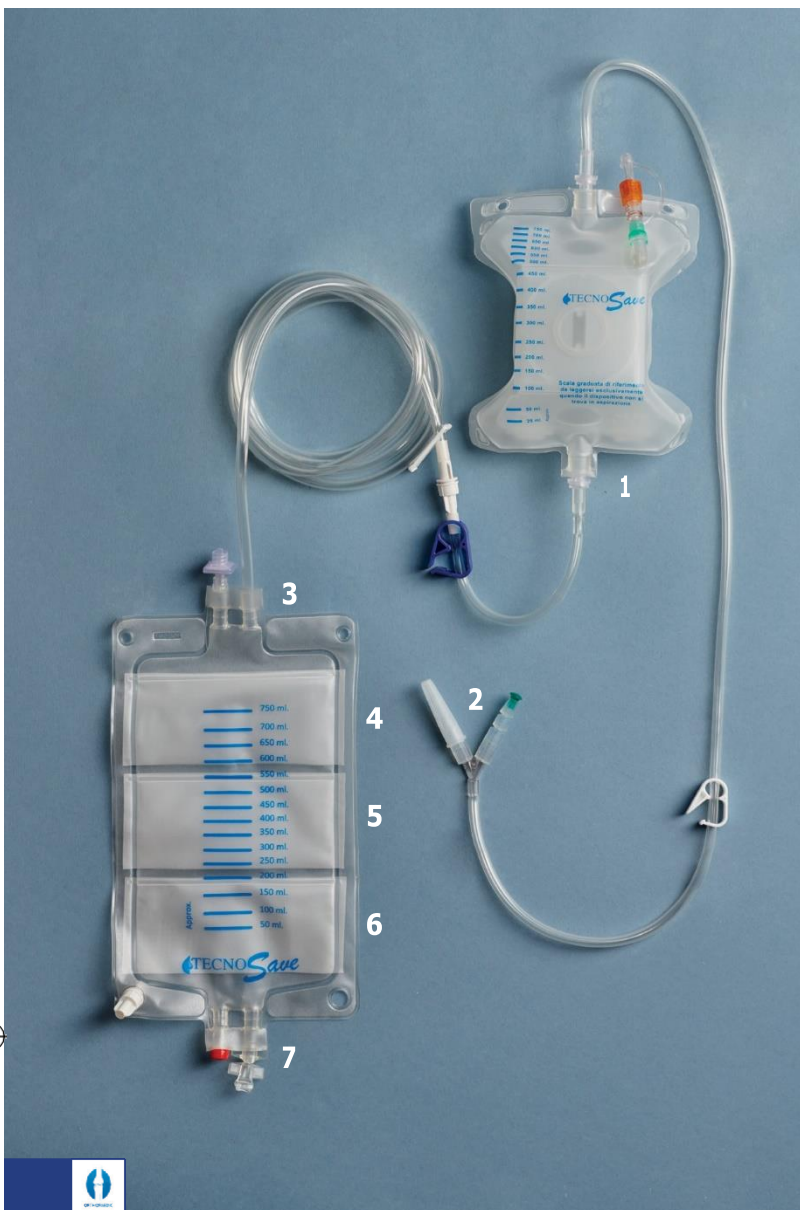

DESCRIPTION

TecnoSave is a complete disposable closed-circuit system for post-operative surgical drainage, collection and recovery of autologous blood for immediate autotransfusion. TecnoSave has an starting pressure of 140 mmHg (18.66 KPa) and an average working pressure of about 40 mmHg (5.33 KPa). The device is provided pre-charged, and due to its design, whenever blood is transferred to the autotransfusion unit for transfusion, it is automatically recharged. The autotransfusion unit consists of a bag with three different levels of filtration : 120 micron macrofilter, 80 micron intermediate filter and a 40 micron micro-aggregate filter.

- 1 - SELF-ACTIVATING COLLECTION UNIT
750ML CAPACITY
- 2 - DOUBLE MULTI-CALIBRE CONNECTOR
- 3 - COLLECTION BAG WITH CASCADE FILTER
- 4 - 120 MICRON FILTERING NET
- 5 - 80 MICRON FILTERING NET
- 6- 40 MICRON FILTERING NET
- 7 - REINFUSION CONNECTION

COD.

DESCRIPTION

QTY/PKG.

TEC859

Orthopaedic post-operative autotransfusion device

10 pc

CE 0051



Tecnoline

Tecnoline S.p.A.
Registered Office Via Gelatti 11/A - 41033 Concordia sulla Secchia - Modena - Italy
Headquarters Via Gelatti 11/A and 18 - 41033 Concordia sulla Secchia - Modena - Italy
T 0535 40312 **M** info@tecnolinemed.com