


SCAN ME

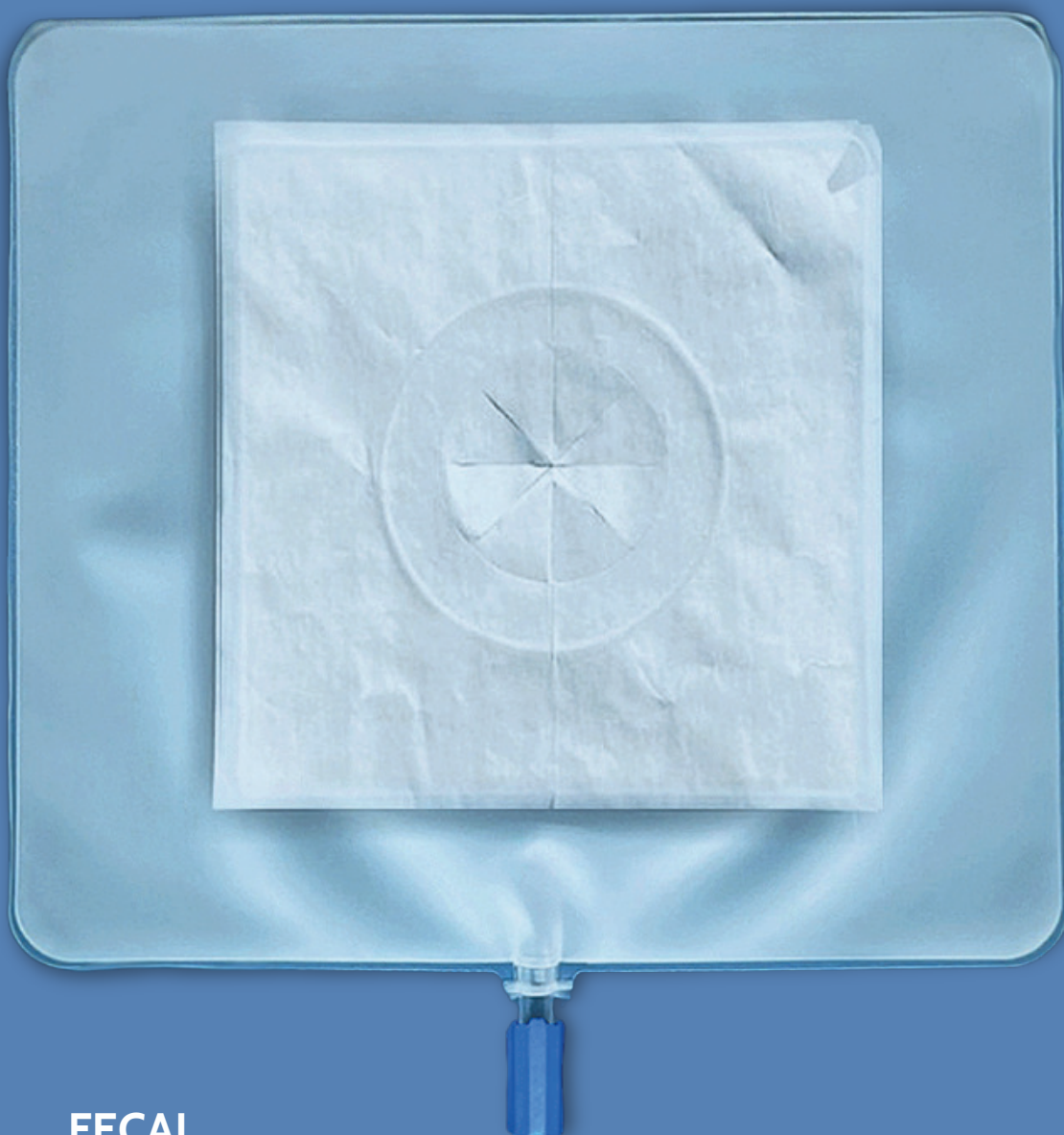

## FECAL CONTAINMENT DEVICE

IGIENITEC is a non-invasive medical device designed for the containment of fecal matter.

It has been developed to prevent environmental and skin contamination caused by fecal soiling.

## FIELDS OF APPLICATION

### FIELDS OF APPLICATION

- **HOSPITAL USE:**

Post-operative management to ensure maximum cleanliness.

- **HOME USE:**

For bedridden individuals or those with limited physical and mental self-sufficiency and/or urinary and fecal incontinence.

- **DAILY MANAGEMENT:**

For active and self-sufficient individuals affected by intestinal disorders with diarrhea and insufficient sphincter control.

IGIENITEC consists of a biocompatible, flexible medical adhesive with a central opening.

The central adhesive part adheres to the perianal area, while the other end is connected to a collection bag.

The distal part of the device includes an access port for inserting a catheter to perform enemas with maximum hygiene and environmental safety.

This access is sealed with a vented cap, allowing the release of flatulence, filtering unpleasant odors, and preventing the bag from inflating due to trapped gas.

LATEX-FREE  
PHTHALATE-FREE

## ORDERING GUIDE

| CODE        | DESCRIPTION | QUANTITY |
|-------------|-------------|----------|
| TEC353608BL | IGIENITEC   | 25 PZ    |



**Tecnoline S.p.A.**  
Registered Office Via Gelatti 11/A - 41033 Concordia sulla Secchia - Modena - Italy  
Headquarters Via Gelatti 11/A and 18 - 41033 Concordia sulla Secchia - Modena - Italy  
T (+39) 0535 40312 M info@tecnolinemed.com



Rev.5 17/07/2024