




Tecnoline

Medical & Surgical Solutions

TECNOSAVE

Post-operative blood
recovery system

SCAN ME

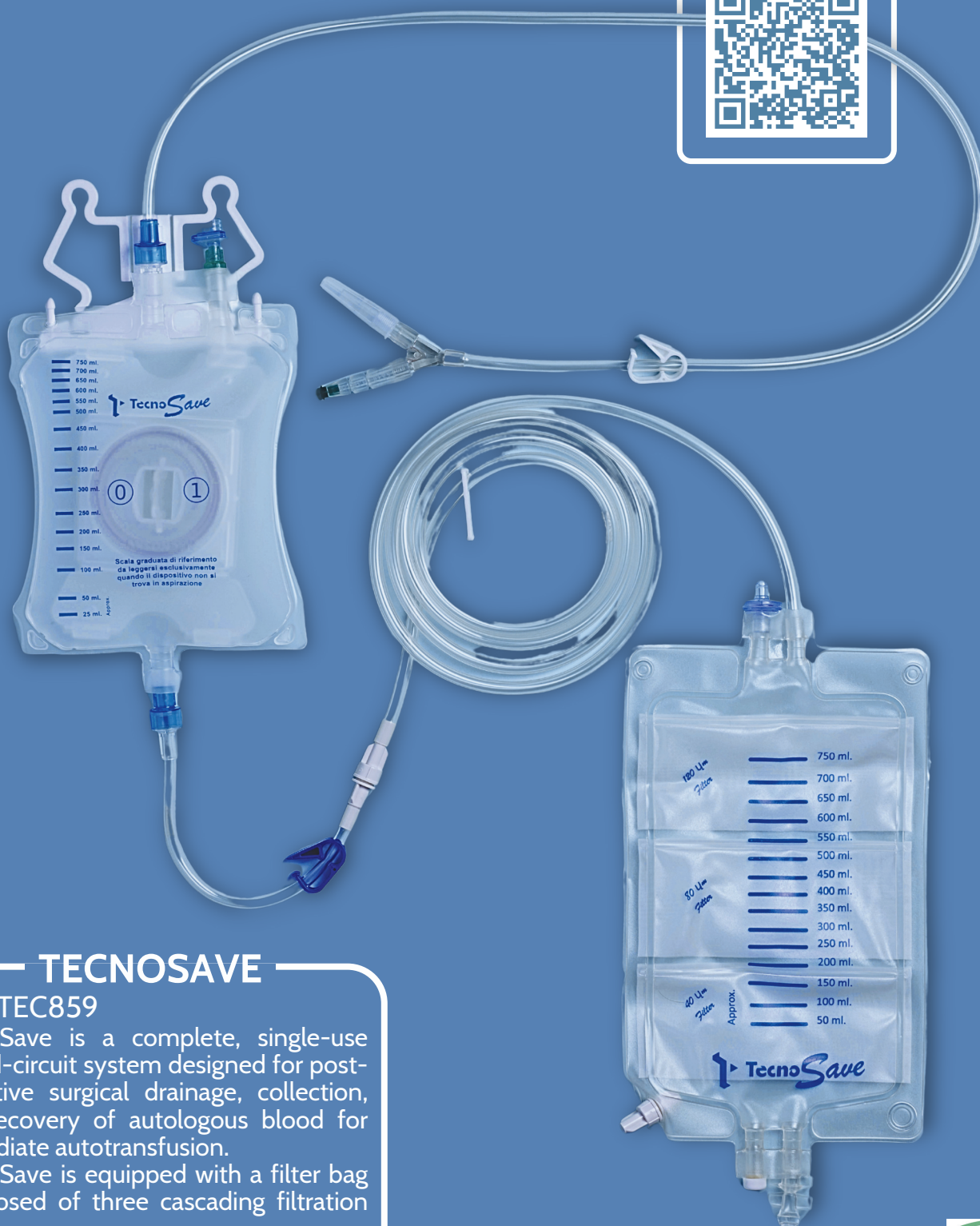

TECNOSAVE

Cod. TEC859

TecnoSave is a complete, single-use closed-circuit system designed for post-operative surgical drainage, collection, and recovery of autologous blood for immediate autotransfusion.

TecnoSave is equipped with a filter bag composed of three cascading filtration levels:

- Macro filter: 120 microns
- Intermediate filter: 60 microns
- Microaggregate filter: 40 microns



Rev.5 17/07/2024

TECNOSAVE SYSTEM

TecnoSave is a complete, single-use closed-circuit system for post-operative surgical drainage, collection, and recovery of autologous blood for immediate autotransfusion.

The device has an initial pressure of 140 mmHg and an average working pressure of approximately 40 mmHg.


TecnoSave is pre-charged and, thanks to its structure, after transferring the blood to the filter bag, it can reactivate itself automatically.

The filter bag includes three cascading filtration levels:

- Macro filter: 120 microns
- Intermediate filter: 60 microns
- Microaggregate filter: 40 microns

TECNOSAVE SYSTEM COMPONENTS

1. Dual multi-caliber connector
2. Self-reactivating 750 ml collection unit
3. Collection bag with cascading filtration
4. 120-micron filter mesh
5. 80-micron filter mesh
6. 40-micron filter mesh
7. Reinfusion connection

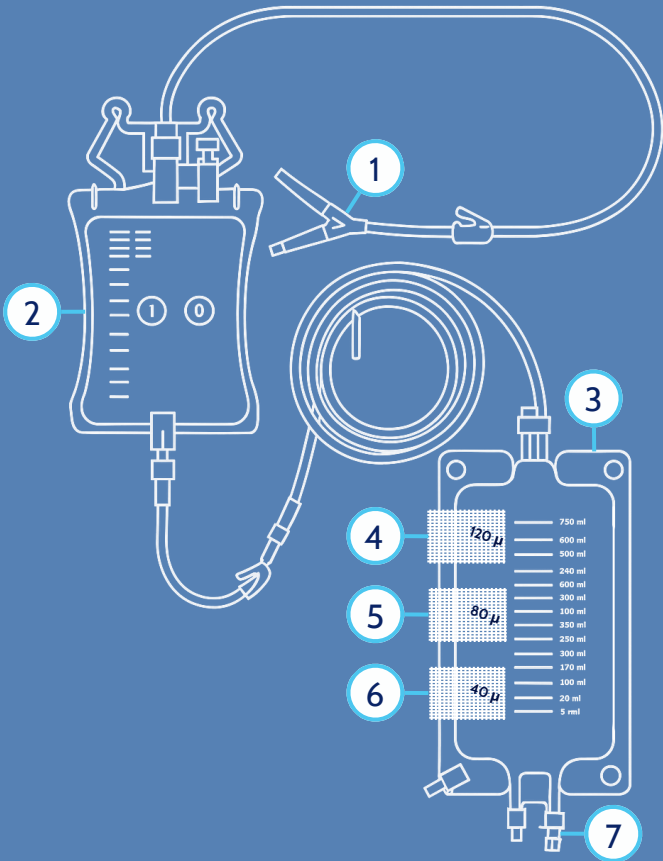

ORDERING GUIDE

CODE	DESCRIPTION	QUANTITY
TEC859	TECNOSAVE	10 PZ



Tecnoline S.p.A.
Registered Office Via Gelatti 11/A - 41033 Concordia sulla Secchia - Modena - Italy
Headquarters Via Gelatti 11/A and 18 - 41033 Concordia sulla Secchia - Modena- Italy
T (+39) 0535 40312 M info@tecnolinemed.com

